





**MUREX**  
DIAGNOSTIC PRODUCTS SPECIALISTS



# URIT UC-58

URINALYSIS CATEGORY

**RELIABLE AND COMPACT DEVICE**  
Semi-Automatic Urine Analyzer

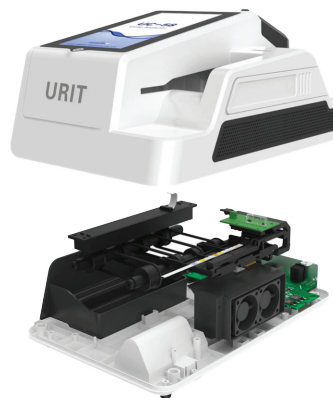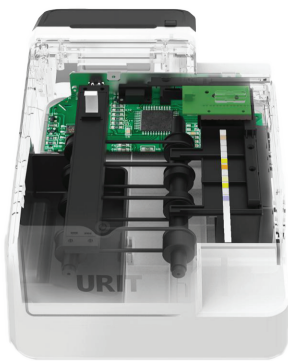

- Up to 550 T/H
- 7-inch extra 50% resolution sensitive touch screen
- Available from 1 to 14 parameters including CR, MA, CA, VitC
- Real strip image display
- Temperature control system integrated
- Lithium battery enhanced

**MUREX** Diagnostic Products Specialists



ISO 9001:2015  
**CERTIFIED**

08 Bonifacio Street, Poblacion East,  
San Nicolas 2447, Pangasinan &  
119 14th Ave., Cubao, Quezon City, Metro

+639 189647158  
+075-600-8159  
+02-8-962-8743

@murexdiagnostics.ph  
www.murexdiagnostics.com  
info@murexdiagnostics.com



SCAN ME



SCAN ME

# UC-58 SPECIFICATIONS

<b>Parameters</b>	18 Parameters Leukocyte, Ketone, Nitrite, Urobilinogen, Bilirubin, Protein, Glucose, Specific Gravity, Blood, pH, Ascorbic Acid, Microalbumin, Creatinine, Calcium, ACR, PCR, Color, Turbidity
<b>Measurement Principle</b>	Reflectance photometry
<b>Reaction Time</b>	45 seconds
<b>Throughput</b>	Up to 550 T/H
<b>Display</b>	7-inch high resolution (1024*600) capacitive color touch screen
<b>Memory Capacity</b>	20,000 results
<b>Strip</b>	HC series strips
<b>Interface</b>	Serial port (Rs232), USB (barcode reader), wifi, bluetooth optional
<b>Work Condition</b>	Environment temperature: 5~30°C, RH≤80%, lamination≤2000lx
<b>Power Supply</b>	Instrument: 12V DC; Adapter: 100V–240V, 50/60Hz
<b>Dimension &amp; Weight</b>	308 mm (L) × 210 mm (W) × 151 mm (H), 3 kg



## Flexible Test Mode

- Continuous measuring 550 T/H
- Strip automatic unloading



## 14 Items Strips

- 14 test items including CR, MA, CA, VitC
- Albumin-to-Creatinine Ratio (ACR) helps to screen early kidney disease
- Color and turbidity results integrated



## Complete QC System

- URIT QC solution provides all test items, including VitC



## Advanced CIS System

- High accuracy and high efficiency
- 5 wavelengths guarantee precise and confident results
- Real image display



## Powerful

- DC 12V & Lithium battery integrated
- More than 600 tests available without power supply



## Easy to Use

- Touchscreen with linux system
- Analyzer automatically detects strip existence, type and starts tests



## Constant Temperature Control System

- Avoid temperature interference
- Allow accurate testing in minimum to 5°C room temperature
- Offer the best reaction environment

MUREX Diagnostic Products Specialists



ISO 9001:2015  
CERTIFIED

08 Bonifacio Street, Poblacion East,  
San Nicolas 2447, Pangasinan &  
119 14th Ave., Cubao, Quezon City, Metro

+639 189647158  
+075-600-8159  
+02-8-962-8743

@murexdiagnostics.ph  
www.murexdiagnostics.com  
info@murexdiagnostics.com



SCAN ME



SCAN ME