

# PORTABLE COAGULATION ANALYZER MC100

COAGULATION CATEGORY

## FAST, SIMPLE, ACCURATE

Optical POCT Coagulation Detection



**SINGLE PERSON TEST DESIGN:** Pre-packaged lyophilized reagents, multiple combo tests of lyophilized reagents save the cost



**SMALL SIZE:** Easy to carry and transport, suitable for multiple scenarios



**VISUALIZATION RESULTS:** Reaction curve & results displayed, more intuitive for clinical judgment



**MULTIPLE METHODOLOGIES:** Optical coagulation method, immunoturbidimetry, chromophoric substrate



**MUREX**  
DIAGNOSTIC PRODUCTS SPECIALISTS

8 Bonifacio Street., Poblacion East, San Nicolas 2447, Pangasinan, PHILIPPINES  
119 14th Ave., Cubao, Quezon City, Metro [info@murexdiagnostics.com](mailto:info@murexdiagnostics.com)  
+639 189647158 | +075-600-8159 | +075-632-4435 | +02-8-962-8743  
[www.murexdiagnostics.com](https://www.murexdiagnostics.com)  
[@murexdiagnostics.ph](https://www.facebook.com/murexdiagnostics.ph)

## Technical specifications

Test parameters	PT, APTT, FIB, TT, DD, AT, FDP, ACT, Anti-Xa, ect.
Sample type	Whole blood, Plasma
Reagent type	Lyophilized reagent
Test time	≤12min
Reagent recognition	Automatic recognition of external QR code on reagent disk
Printer	Built-in thermal printer, also available to external printer ( with wireless printer drive)
Operating system	Android 5-inch( Version 5.1), distinguishability 720×1280, support LIS/HIS system
Size	18*27*15.5cm
Weight	≤3KG

## Operation Process



Put in the disk

1




Add sample

2



Click 'start'

3

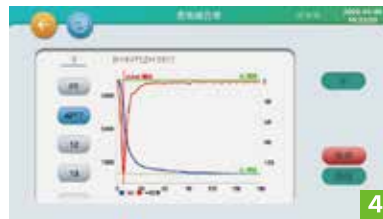

Read results

4

## Applicable to



ICU



Maternity & children



Dentistry



Ophthalmology



Clinical departments



Orthopedics



Field first aid



Hospital specialties



Ambulance



Zhejiang Pushkang Biotechnology Co., Ltd.

Science & Technology Park, No.398 Mahuan Road, Binhai New Area, Shaoxing, Zhejiang.

Hotline: 400-003-9660

Email: biz@pushkang.cn

Web: www.pushkang.com.cn

